


# PRODUCT DATA SHEET

## Flexima® 3S Colostomy Bags Two-piece stoma bag systems

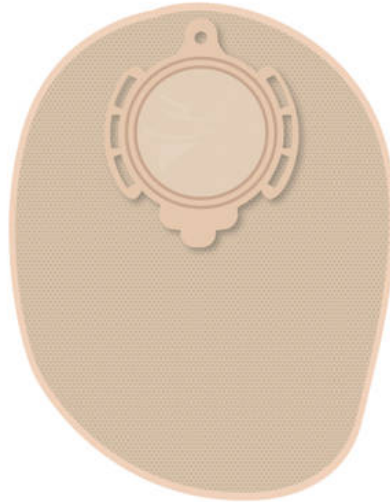

Flexima® 3S is a 2-piece range with a guided mechanical coupling system for an easy positioning of the stoma bag onto the wafer. Its ergonomically shaped colostomy bags are available in two different sizes, midi and maxi.

### Colors and capacity (volume as per ISO 8670-2):

- **Midi:**
  - Ø45 mm: capacity of +/- 520 ml, beige
  - Ø55 mm: capacity of +/- 500 ml, beige
  - Ø65 mm: capacity of +/- 465 ml, beige
- **Maxi:**
  - Ø45 mm: capacity of +/- 665 ml, beige and transparent
  - Ø55 mm: capacity of +/- 655 ml, beige and transparent
  - Ø65 mm: capacity of +/- 630 ml, beige and transparent
  - Ø80 mm: capacity of +/- 510 ml, beige and transparent

### High performance protected filter (HP):

- Improved deodorization capacity through its high-tech breathable membrane
- Leak-proof and permeable to gases

### Material:

- Soft non-woven and odor-proof plastic film
- Free of latex (PVC) and phthalates (DEHP)

**Others:**

- To be used with the Flexima® 3S wafer corresponding to the chosen diameter of the bag
- Box of 30 bags

**Advantages**

Flexima® 3S stand for Simple, Soft and Secure:

- **Simple** coupling system that intuitively guides you through the attachment process so that only a minimum of patient education is required. Depending on which of the three slots you choose a multiple bag orientation is possible to adapt the bag positioning to your body contours.
- **Soft** coupling system with a floating flange that requires a minimum of pressure on the abdomen to prevent discomfort during application. The coupling system remains flexible when the wafer and bag are coupled together.
- **Secure** coupling system for three reasons: Enhanced leak-proof thanks to the close fitting of the bag on the wafer, increased security through the click when wafer and bag are well locked together and finally even more fixation with the possibility to use a belt attached to one of the three hooks on the bag.

Article Description	Article Number
FLEXIMA 3S COLO MIDI BEIGE 45MM	931045N
FLEXIMA 3S COLO MIDI BEIGE 55MM	931055N
FLEXIMA 3S COLO MIDI BEIGE 65MM	931065N
FLEXIMA 3S COLO MAXI BEIGE 45MM	931445N
FLEXIMA 3S COLO MAXI BEIGE 55MM	931455N
FLEXIMA 3S COLO MAXI BEIGE65MM	931465N
FLEXIMA 3S COLO MAXI BEIGE 80MM	931480N
FLEXIMA 3S COLO MAXI TRP 45MM	931545N
FLEXIMA 3S COLO MAXI TRP 55MM	931555N
FLEXIMA 3S COLO MAXI TRP 65MM	931565N
FLEXIMA 3S COLO MAXI TRP 80MM	931580N